

# Determination of Phenoxy-acid Herbicides in Various Matrices

Thermo Fisher Scientific Inc., Austin, TX, USA

## Introduction

Phenoxy-acid herbicides are often used mixed with fertilizers and applied to lawns and gardens to remove dandelions, thistles, and ragweed. The herbicides are toxic to these broadleaf plants, which are absorbed through the plant surface and transported throughout the circulatory system to all parts of the plant. The plant dies due to root starvation caused by abnormal cell growth that blocks the movement of liquids and nutrients. The herbicides will exist in the garden for up to three weeks exposing other animals to the effects of this poison. Also, the runoff of the lawns and gardens can pose a threat to the animals of nearby rivers and lakes. Since the long-term effects of these compounds are not completely known, monitoring for the presence of phenoxy acid herbicides is important.<sup>1</sup>

The Thermo Scientific PolarisQ external source ion trap mass spectrometer operated in EI-MS/MS will be used to demonstrate a sensitive and robust methodology for the analysis of phenoxy acid herbicides in various matrices. This is accomplished by coupling the splitless injection technique with MS/MS. This ensures that both sensitivity and confirmation can be achieved in pure standards as well as matrix.

## Instrument Conditions

### PolarisQ Ion Trap

Ion Source temperature: 250 °C  
Ionization Mode: +EI; 70 eV  
AGC: 50  
Injection Waveforms: On (default)  
Buffer Gas Flow: 2.0 mL/min  
MS/MS Parameters: See Table 1 on page 3.

### TRACE GC Ultra Gas Chromatograph

Column: Restek Rtx™-5; 0.18 mm ID x 40 m; 0.2 µ film thickness  
Oven: 100 °C (1 min); 10 °C/min  
300 °C (4 min); 25 min run time  
Split/Splitless injector: Splitless mode  
Injector temperature: 250 °C  
Column flow: 1.5 mL/min  
Splitless time: 1.0 min  
Injection port liner: 5 mm Focus Liner; splitless; packed with glass wool

### Autosampler

Injection Volume: 2 microliters

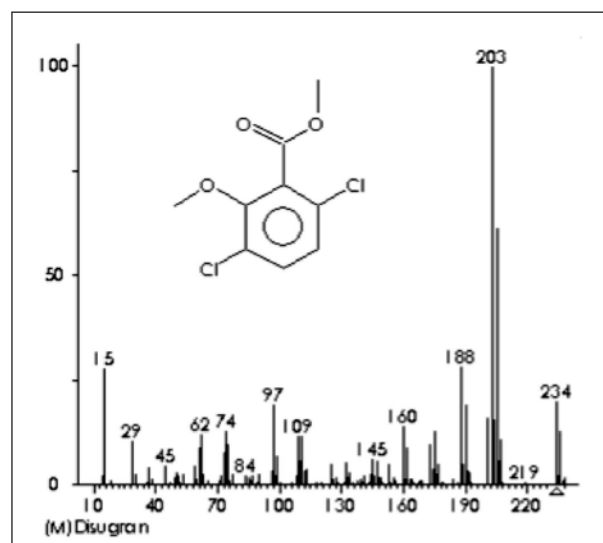


Figure 1: Full scan EI spectrum of dicamba methyl ester from NIST library.

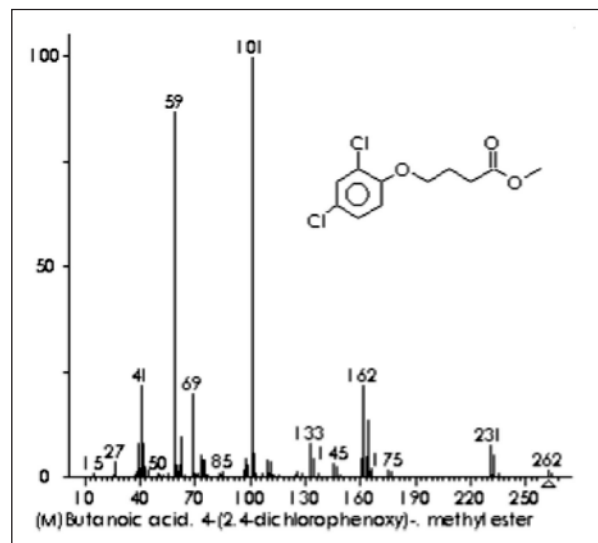


Figure 2: Full scan EI spectrum of 2,4-DB methyl ester from NIST library.

## Key Words

- PolarisQ
- Environmental
- GC-MS/MS
- Herbicides
- Ion Trap
- Water Quality

## Objective

The purpose of this study was to develop a robust and sensitive method for the determination of methylated phenoxy acid herbicides by GC-MS/MS. Utilizing USEPA Method 8151 as the starting point, a calibration curve containing herbicides was purchased from AccuStandard. The method was evaluated for three matrices, soil, pond water, and lawn clippings. Each matrix was injected multiple times to validate detectivity and reproducibility.

## Sample Preparation

500 mL of pond water was extracted twice with 60 mL of methylene chloride. 50 grams of soil was mixed with 50 grams of sodium sulfate and sonicated for 15 minutes with 100 mL methylene chloride. 10 grams of lawn clippings were mixed with 15 grams sodium sulfate and sonicated with 150 mL methylene chloride for 15 minutes. The methylene chloride extracts were dried with sodium sulfate and the resulting extracts were evaporated to dryness with nitrogen. The prepared extracts were reconstituted with 2000  $\mu$ L of isooctane. 900  $\mu$ L of the isooctane was transferred into two vials. 100  $\mu$ L of isooctane was added to the blank vial and 100  $\mu$ L of the methylated standard "Level 6" was added to the "spike" vial.

## Results and Discussion

Electron impact was chosen because of the consistent response of all compounds. The MS/MS parameters for each compound are listed in Table 1 and the fragmentation pathways for two of the compounds are shown in Figures 3 and 4. In each case the base peak was chosen to be the precursor ion. See Figures 1 and 2 on page 1 for dicamba methyl ester and 2,4-DB methyl ester.

A calibration curve was prepared by injecting each standard in triplicate. See Figures 5 and 7 for the mass chromatograms and Figures 6 and 8 for linearity graphs of dicamba methyl ester and 2,4-DB methyl ester respectively. The chromatograms in Figure 9 show the total ion chromatograms from the full scan analysis of the various extracts before spiking. Table 1 summarizes the calibration range and correlation coefficients for all compounds. Each of the three matrices was injected 12 times for reproducibility. External standardization was used for calibration purposes to emulate a worst-case scenario from the standpoint of system robustness, since it will not compensate for the loss of signal due to ion source or any other contamination. By injecting spiked matrix samples repetitively and plotting the resulting areas over time one can surmise how well a particular GC-MS system will perform a given method. In this case we have chosen to extract and inject soil, lawn clippings, and pond water without regard to any sample clean-up. The confirmation of peaks was performed by forward library searching against a spectrum obtained from a standard (see Table 2).

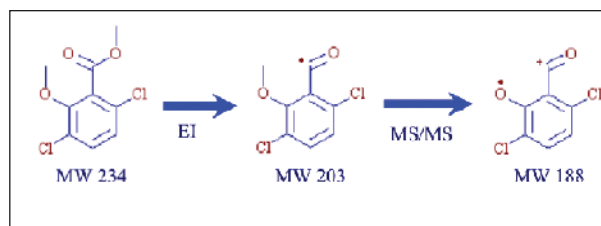


Figure 3: The EI and subsequent MS/MS fragmentation of dicamba methyl ester.

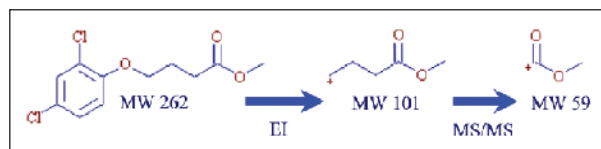


Figure 4: The EI and subsequent MS/MS fragmentation of 2,4-DB methyl ester.

COMPOUND	RETENTION TIME	PRECURSOR ION	QUAN.ION	CONC.RANGE (pg/ $\mu$ L)	CORRELATION COEFF.
Dicamba	10.15	203	188	10-350	0.9974
MCPP	10.49	228	169	2000-100,000	0.9935
MCPA	10.70	214	182	2000-100,000	0.9948
Dichlorprop	11.21	162	126	20-700	0.9963
2,4-D	11.45	199	125	20-700	0.9950
2,4,5-TP	12.83	196	132	5-175	0.9977
2,4,5-T	13.14	232	159	5-175	0.9990
2,4-DB	13.82	101	59	20-700	0.9963
Dinoseb	19.90	225	195	5-175	0.9937

Table 1: MS/MS Parameters and Calibration Results

COMPOUND	SPIKE			GRASS			POND WATER			SOIL		
	pg/ $\mu$ L	FIT STD.		pg/ $\mu$ L	%RSD	FIT STD.	pg/ $\mu$ L	%RSD	FIT STD.	pg/ $\mu$ L	%RSD	FIT STD.
Dicamba	35	987		31	3.2%	911	31	2.2%	922	32	4.3%	916
MCPP	6000	937		5034	4.6%	880	5367	3.5%	887	5292	5.1%	877
MCPA	6000	973		4809	2.6%	917	5131	3.5%	920	5024	5.4%	917
Dichlorprop	70	995		66	3.1%	902	64	3.5%	925	65	4.5%	914
2,4-D	70	987		69	2.5%	897	65	3.2%	896	66	4.4%	896
2,4,5-TP	17	992		17	3.5%	750	17	4.0%	853	17	6.0%	777
2,4,5-T	17	995		18	4.6%	690	18	4.9%	777	18	6.8%	695
2,4-DB	70	1000		68	3.7%	892	69	3.9%	918	69	5.1%	911
Dinoseb	17	979		28	8.0%	820	24	8.3%	921	27	13.7%	899

Table 2: Matrix Spike Recovery Data

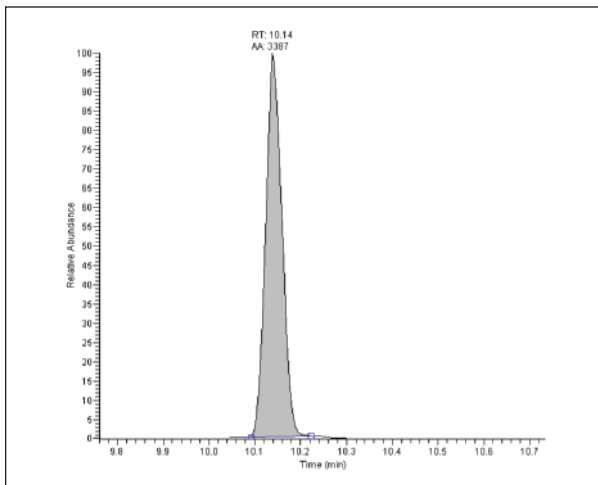


Figure 5: Mass chromatogram of MS/MS product ion for dicamba methyl ester at 10 pg/ $\mu$ L

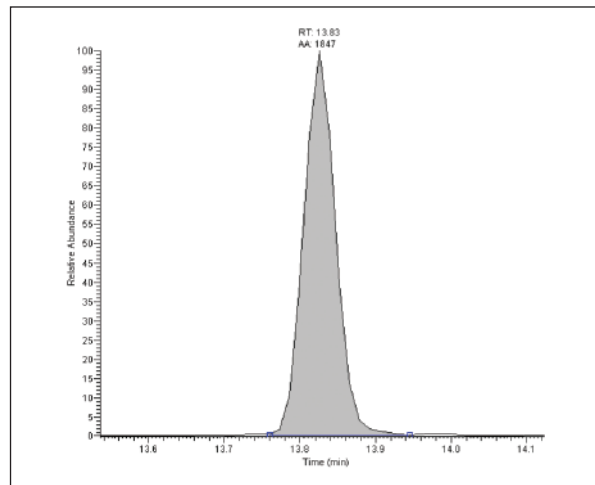


Figure 7: Mass chromatogram of MS/MS product ion for 2,4-DB methyl ester at 20 pg/ $\mu$ L

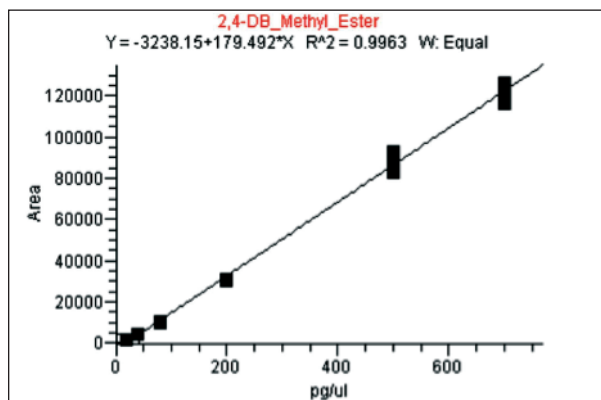


Figure 6: Calibration curve for dicamba methyl ester from 10 pg/ $\mu$ L to 350 pg/ $\mu$ L

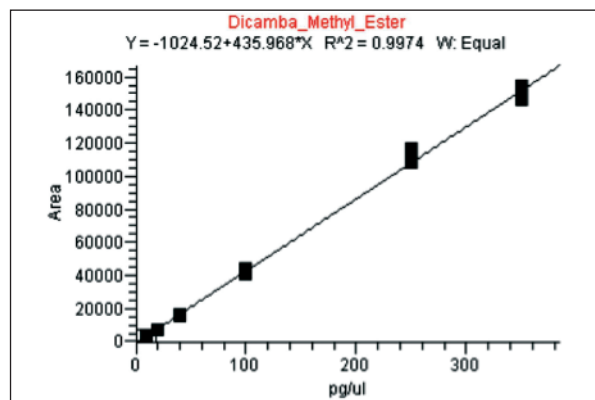


Figure 8: Calibration curve for 2,4-DB methyl ester from 20 pg/ $\mu$ L to 700 pg/ $\mu$ L

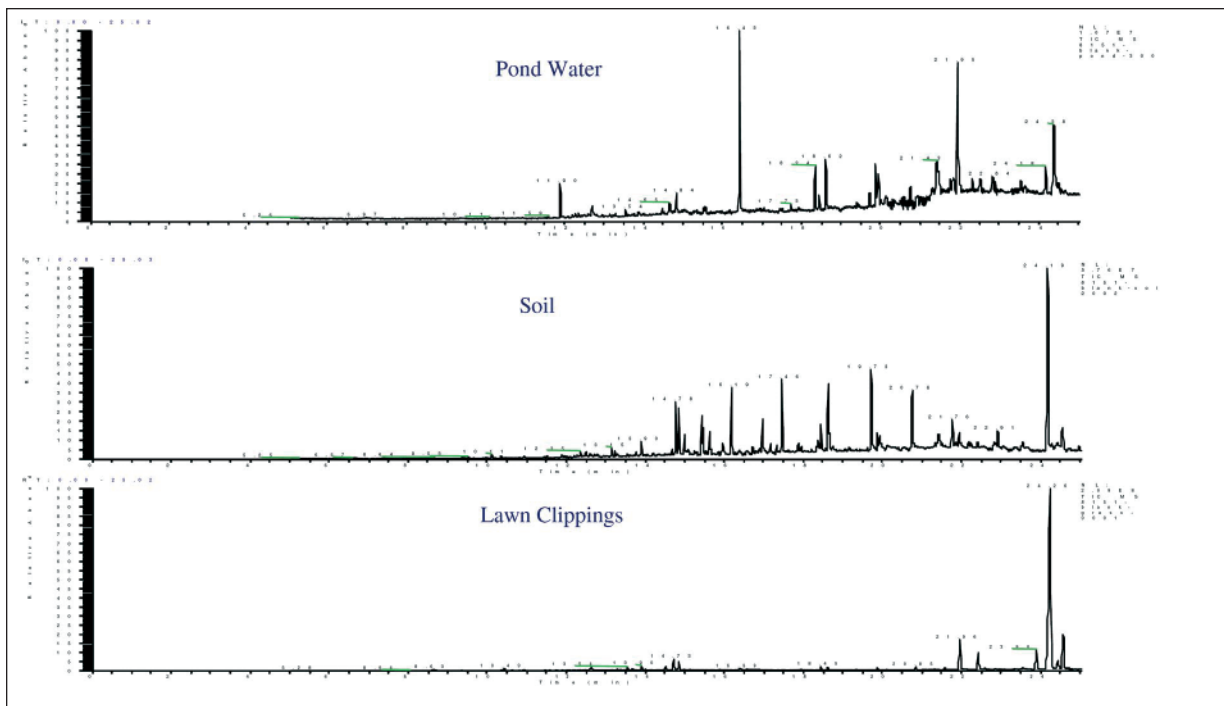


Figure 9: Full scan chromatograms of the pond water, soil, and lawn clippings extracts before spiking.

## Conclusion

Although these compounds are normally analyzed by GC-ECD the use of GC-MS/MS eliminates the additional confirmation requirement of a secondary column and/or additional GC-MS analysis<sup>2</sup>. From the results in Table 1 the method is linear and quantitative. This was all performed without significant sample cleanup or any maintenance of the system during the entire method development and sample analysis. The use of MS/MS also provides a level of confirmation not available with conventional GC detectors. From the forward spectral library matching results (Fit - where a value of 1000 is a perfect match) in Table 2 it is clear that the methodology (MS/MS) is not only robust and quantitative but it is sensitive and provides confirmation in a single injection without the need for extensive sample preparation.

## References

1. <http://www.mb.ec.gc.ca/pollution/pesticides/ec00s10en.html>
2. Method 8151A, Chlorinated Herbicides by GC Using Methylation or Pentfluorobenzoylation Derivatization US EPA, Revision 1, December 1994

## Acknowledgement

Authors: John D. Ragsdale III and Meredith Conoley

In addition to these offices, Thermo Fisher Scientific maintains a network of representative organizations throughout the world.

**Australia**  
+61 2 8844 9500

**Austria**  
+43 1 333 50340

**Belgium**  
+32 2 482 30 30

**Canada**  
+1 800 532 4752

**China**  
+86 10 5850 3588

**Denmark**  
+45 70 23 62 60

**France**  
+33 1 60 92 48 00

**Germany**  
+49 6103 408 1014

**India**  
+91 22 6742 9434

**Italy**  
+39 02 950 591

**Japan**  
+81 45 453 9100

**Latin America**  
+1 608 276 5659

**Netherlands**  
+31 76 587 98 88

**South Africa**  
+27 11 570 1840

**Spain**  
+34 91 657 4930

**Sweden/Norway/Finland**  
+46 8 556 468 00

**Switzerland**  
+41 61 48784 00

**UK**  
+44 1442 233555

**USA**  
+1 800 532 4752

[www.thermo.com](http://www.thermo.com)



Thermo Fisher Scientific Inc., Austin, Texas, USA is ISO Certified.

©2007 Thermo Fisher Scientific Inc. All rights reserved. Siltek and Rtx are registered trademarks of Restek Corporation. All other trademarks are the property of Thermo Fisher Scientific and its subsidiaries.

Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.

AN10041\_E 08/07C